



“Mobility, Comfort and Style”



X Series

www.aktiv-wheelchairs.com



The Aktiv Series range of manual wheelchairs are designed to provide **comfort, safety** and a **stylish** and contemporary look with a practical balance of durability and functionality in mind.

DESIGN

We've endeavoured to create a range that gives a choice of quality wheelchairs to suit a wide variety of needs from basic occasional use to the more specific requirements of a long term user or person with more demanding physical needs. We have listened intently to the feedback of users, carers and professionals to help us to produce a range which is **reliable, high quality, multifunctional** and **economical**.

HEALTH SERVICE FRIENDLY

We have considered the needs of the health services who may provide and reuse these wheelchairs. Therefore we have designed our chairs so that most **parts are interchangeable** between models e.g. a leg rest or armrest assembly for an X3 with also fit an X4, X5, X6 and our accessories include trays, oxygen carriers, drip stands and amputee supports also fit easily.

QUALITY

We have used **quality fittings** for moving parts which are subject to high wear e.g. the clamps we use on our armrest fittings are chrome rather than plastic and where we do use plastic it is in low stress areas and only high quality virgin plastic which is **extremely strong** is used. We fit all our chairs with PU tyres which are **puncture proof** as standard to avoid costly and inconvenient punctures however pneumatic tyres are available on request.

X1 Series - Basic Lightweight Wheelchair



The **Aktiv X1** is a fantastic value and durable wheelchair for occasional use. This wheelchair has a strong steel folding frame for easy storage and transportation. The Armrests are fixed position as standard and have comfortable padded armrest pads. The backrest folds down and the leg rests and rear wheels are removable by easy quick release mechanisms, therefore the chair can be quickly and easily disassembled for easy storage in a car boot. This chair has the added benefit of a dual axle meaning that the chair can be easily and economically converted from a self-propel to an attendant transit wheeled model.

- Colour: Metallic Sky Blue
- Seat Width: 16" (40cm) , 18" (45cm), 20" (50cm)
- Seat Depth: 16" - 18" (40-45cm)
- Seat Height: 20" (50cm)
- Backrest Height: 17" (43cm)
- Total Weight: 9.2kg
- Transport Weight: 8kg

FEATURES

- Steel folding frame
- Fold down Backrest
- Available in Transit (12" Rear wheel) or Self-Propel Models (24" rear wheel)
- Dual axle allows easy conversion between transit & self propel models
- Available as fixed non-removable armrest or flip back removable armrests with padded arm pad
- Steel side panel clothes guard
- Swing away leg rests with black plastic foot plate & heel loops
- Black nylon seat and back upholstery
- 7" grey PU front castors
- Choice of press button or buckle style Safety belt

Armrest Options page 20.

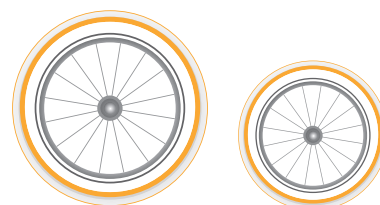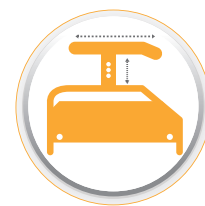

Quick release wheels: 24" - 12"



Dismantle for transport.



Easy fold away.



Swing away detachable leg rest.



**Max User Weight:
121kg/19st.**



Dual axle allows easy interchange of TR/SP wheels.



Crash Tested



Puncture Proof Solid Tyres

Product complies with EN12183 and Crash Test ISO7176-19:2001.

X2 Series - Lightweight Aluminium Wheelchair

The **Aktiv X2** is a lightweight durable wheelchair for occasional use. This wheelchair has a lightweight aluminium folding frame and is easy to lift, store and transport. The **armrests flip back and can be removed** easily as standard and have comfortable polyurethane armrest pads – Interchangeable armrest styles available. The backrest folds down and the leg rests and rear wheels are removable by easy quick release mechanisms; therefore the chair can be quickly and easily disassembled for easy storage in a car boot. This chair has the added benefit of a **dual axle** meaning that the chair can be easily and economically converted from a self-propel to an attendant transit wheeled model.



- Colour: Metallic Sky Blue
- Seat Width: 16" (40cm) , 18" (45cm), 20" (50cm)
- Seat Depth: 16" - 18" (40-45cm)
- Seat Height: 20" (50cm)
- Backrest Height: 17" (43cm)
- Total Weight: 13.8kgs
- Transport Weight: 7.2kgs

Product complies with EN12183 and Crash Test ISO7176 - 19:2001

FEATURES

- Lightweight Aluminium wheelchair
- Available in Transit (12" Rear wheel) or Self-Propel Models (24" rear wheel)
- Dual axle allows easy conversion between transit & self propel models
- Aluminum frame
- Flip back and detachable armrest with PU arm pad
- Plastic side panel clothes guard
- Swing away leg rest with black plastic foot plate & heel loop
- Black nylon seat and back upholstery
- 6" grey PU front castors
- Choice of press button or buckle style Safety belt

Interchangeable Armrest
Options page 20.

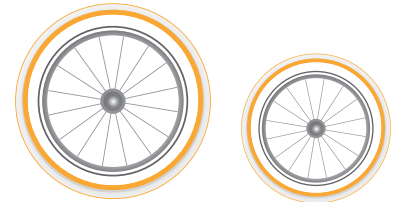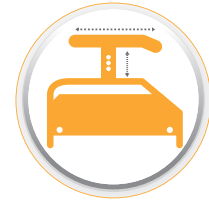

Quick release wheels: 24" - 12"



Easy transportable.



Easy fold away.



Swing away detachable leg rest



**Max User Weight:
121kg/19st.**



Dual axle allows easy interchange of TR/SP wheels.



Crash Tested



Puncture Proof Solid Tyres

X3 Lite - Deluxe Aluminium Wheelchair

The **Aktiv X3 Lite** is a deluxe lightweight durable wheelchair for occasional to full-time use. This wheelchair has a lightweight aluminium folding frame and is easy to lift, store and transport. The armrests flip back and can be removed easily as standard and have comfortable polyurethane armrest pads. The backrest folds down and the leg rests and rear wheels are removable by easy quick release mechanisms; therefore the chair can be quickly and easily disassembled for easy storage in a car boot. This chair **seat height can be adjusted** and can be easily and economically converted from a self-propel to an attendant transit wheeled model if required. The **backrest is tension adjustable** with comfortable padded nylon cover to accommodate and support the users back and a **2" seat cushion**. The chair is fitted with anti-tip wheels for added safety and a stepper to help an attendant to ease the chair up over kerbs.

- Colour: Metallic Sky Blue
- *See Accessories Section for choice of Armrests & Leg rests Available
- Seat Width: 16" (40cm), 18" (45cm), 20" (50cm)
- Seat Depth: 16" - 18" (40-45cm)
- Seat Height: Adjustable 16" (40cm) - 20" (50cm)
- Backrest Height: 17" (43cm)
- Total Weight: 13.8kg
- Transport Weight: 6.8kg



FEATURES

- Folding Aluminum frame
- Fold Down back
- Multiple axle for seat height + depth adjustments
- Flip back and detachable armrest with PU arm pad
- Plastic side panel clothes guard
- Swing away leg rests with black plastic foot plate (Angle adjustable) & heel loop
- Aluminum brake
- 7" grey PU castors
- Anti-tipper
- Choice of press button or buckle style safety belt
- Black nylon seat.
- Black nylon tension back (8 adjustable Velcro straps) with nylon cover
- 2" Seat cushion

Optional:

- Elevating leg rests (see pg20)

Interchangeable Armrest Options page 20.

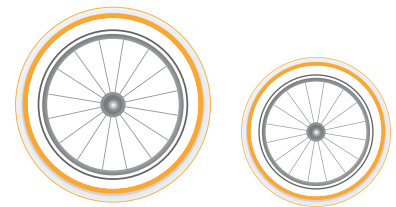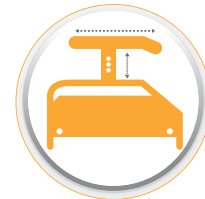

Quick release wheels: 24" - 12"



Easy transportable.



Upholstery with pocket.



**Max User Weight:
130kg/20.5st.**



Crash Tested



Puncture Proof Solid Tyres

Product complies with EN12183 and Crash Test ISO7176-19:2001.

X4 Series - Aluminium Heavy Duty Wheelchair

The **Aktiv X4** is a deluxe **heavy duty** wheelchair for occasional to full-time use for heavier users. This wheelchair has a lightweight aluminium folding frame and is easy to lift, store and transport. The armrests flip back and can be removed easily as standard and have comfortable polyurethane armrest pads. The backrest folds down and the leg rests and rear wheels are removable by easy quick release mechanisms; therefore the chair can be quickly and easily disassembled for easy storage in a car boot. This chair seat height can be adjusted and can be easily and economically converted from a self-propel to an attendant transit wheeled model if required. The **backrest is tension adjustable** with comfortable padded nylon cover to accommodate and support the users back and a **2" seat cushion**. The chair is fitted with **anti-tip wheels** for added safety and a stepper to help an attendant to ease the chair up over kerbs. This chair comes with a **slide on one piece footboard** which fits easily over the standard footplates to give a strong full width platform for heavy users to place their feet comfortably.

- Colour: Metallic Sky Blue
- *See Accessories Section for choice of Armrests & Leg rests Available
- Seat Width: 20" (51cm), 22" (56cm) , 24" (61cm)
- Seat Depth: 17" (43cm)
- Seat Height: Adjustable 16" (41cm) - 20" (51cm)
- Backrest Height: 17" (43cm)
- Total Weight: 12.9kg
- Transport Weight: 8.4kg



FEATURES

- Folding Aluminium Frame
- Fold down back rest
- Multiple axles to allow various seat heights / interchange from transit to self propel wheels
- Flip back and detachable armrest with PU arm pad
- Plastic side panel clothes guard
- Swing away leg rest with black plastic foot plate (angle adjustable) & one piece footplate cover
- Black nylon seat upholstery
- Black nylon tension back with nylon cover (8 adjustable Velcro straps)
- 200x38mm grey PU front castors
- Anti-tipper
- Choice of press button or buckle style safety belt
- Available in Transit or Self Propel models
- 2" Seat cushion

Optional:

- Elevating Leg rests (see pg20)
- Tension adjustable Seat canvas

Interchangeable Armrest Options page 20.

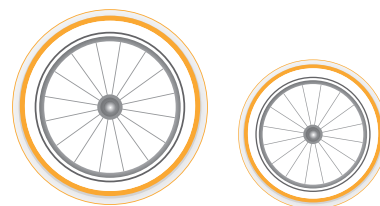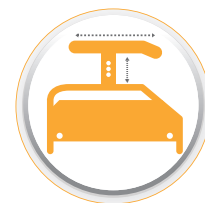

Quick release wheels: 24" - 12"



**Max User Weight:
130kg/20.5st.**



Individual footplates with separate one piece footplate cover.



Removable footplate cover for strength plus easy foot positioning.



Tension adjustable backrest.



Crash Tested

Puncture Proof Solid Tyres

Product complies with EN12183 and Crash Test ISO7176-19:2001.

X5 Series - Deluxe Modular Wheelchair

The **Aktiv X5** is a **modular lightweight** wheelchair suitable for occasional to full-time use and is multi-adjustable. This wheelchair has a lightweight aluminium folding frame and is easy to lift, store and transport. The armrests flip back and can be removed easily as standard and have comfortable polyurethane armrest pads. The backrest folds down and the leg rests and rear wheels are removable by easy quick release mechanisms; therefore the chair can be quickly and easily disassembled for easy storage in a car boot. The modular style frame allows the **seat height, seat depth** and **backrest angle** to be simply and easily adjusted. The chair can be easily and economically converted from a self-propel to an attendant transit wheeled model if required. The backrest is **tension adjustable** with comfortable padded nylon cover to accommodate and support the users back and a **2" seat cushion** is supplied as standard. The chair is fitted with **anti-tip wheels** for added safety and a stepper to help an attendant to ease the chair up over kerbs. **Footplates are height and angle adjustable** as standard.

- Colour: Metallic Sky Blue
- *See Accessories Section for choice of Armrests & Leg rests Available
- Seat Width: 16" (41cm), 18" (46cm), 20" (51cm)
- Seat Depth: 16" (41cm) to 20" (51cm)
- Seat Height: Adjustable 16" (41cm) - 20" (51cm)
- Backrest Height: 17" (43cm)
- Total Weight: 12.9kg
- Transport Weight: 8.4kg



FEATURES

- Folding aluminium frame
- Back height adjustable
- Flip back and detachable armrest with height adjustable PU arm pad
- Plastic side panel clothes
- Swing away leg rest
- Black nylon seat with breathable cover
- Tension adjustable back with wipeable cover
- Back angle adjustable
- 6" grey PU castors
- Available in Transit or Self Propel models
- Choice of press button or buckle style safety belt
- Height adjustable Anti-tipper


Optional:

- Elevating Leg rests (see pg20)

Interchangeable Armrest Options page 20.

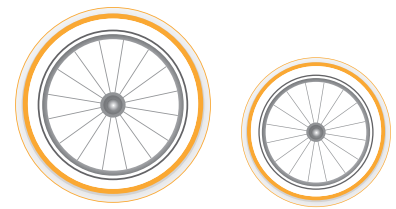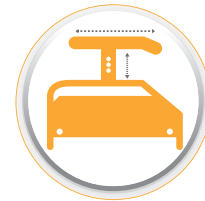

Quick release wheels: 24" - 12"

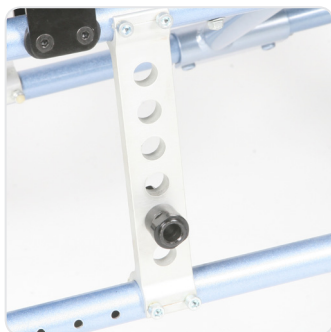

Seat height plus depth adjustable multi axle.



Easy to transport.



Back angle adjustments.



Tension adjustable backrest - 8 straps



**Max User Weight:
130kg/20.5st.**



Crash Tested



Puncture Proof Solid Tyres

Product complies with EN12183 and Crash Test ISO7176-19:2001.

X6 Series - Paediatric Aluminium Wheelchair

The **Aktiv X6** is a **paediatric modular lightweight** wheelchair suitable for occasional to full-time use and is multi-adjustable. This wheelchair has a lightweight Aluminium folding frame and is easy to lift, store and transport. The armrests flip back and can be removed easily as standard and have comfortable polyurethane armrest pads. The leg rests and rear wheels are removable by easy quick release mechanisms; therefore the chair can be quickly and easily disassembled for easy storage in a car boot. The chair can be easily and economically converted from a self-propel to an attendant transit wheeled model if required. The chair is fitted with **anti-tip wheels** for added safety and a stepper to help an attendant to ease the chair up over kerbs. **Footplates are height and angle adjustable as standard** and are supplied with a **one piece footplate** which simply slides onto the standard footplates to give a no-gaps platform for children to place their feet. **Push handles are height adjustable** as standard. The **backrest is tension adjustable** and the upholstery is breathable and padded for extra comfort.

- Colour: Metallic Sky Blue, Princess Pink, Fire Engine Red, Sunshine Yellow
- Seat Width: 12" / 14" / 16"
- Seat Depth: 16" (41cm) to 20" (51cm)
- Seat Height: Adjustable 16" (41cm) - 20" (51cm)
- Backrest Height: 17" (43cm)
- Total Weight: 14.0kg
- Transport Weight: 10.9kg



FEATURES

- Folding aluminium frame
- Push handles height adjustable
- Flip back and detachable armrest with height adjustable PU arm pad
- Plastic side panel clothes
- Swing away leg rest with black plastic foot plate/calf strap/additional one piece footplate cover
- Black nylon seat with breathable cover
- Tension adjustable back with breathable wipeable cover
- Back angle adjustable
- 6" grey PU castors
- Available in Transit or Self Propel models
- Choice of press button or buckle style safety belt
- Height adjustable Anti-tipper

Optional:

- Elevating Leg rests (see pg20)

Interchangeable Armrest Options page 20.

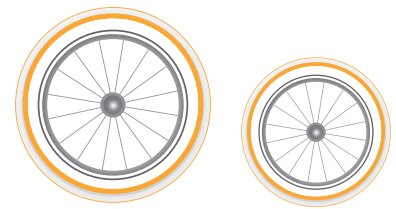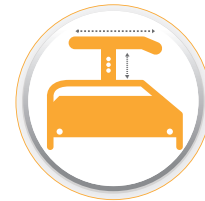

Quick release wheels: 24" - 12"




One piece footplate cover.



Removable cover.



Swing away detachable leg rest: 90° - 170°.

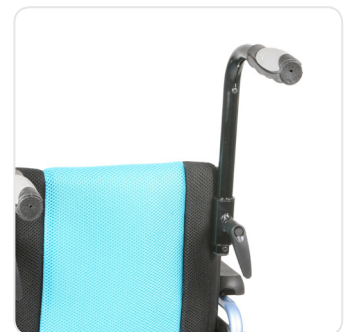

Height adjustable handles as standard.



**Max User Weight:
75kg/12st.**



Crash Tested



Puncture Proof Solid Tyres

Product complies with EN12183 and Crash Test ISO7176-19:2001.



X7 Series - Tilt and Recline Wheelchair

The Aktiv X7 is a **Multi-positional modular** wheelchair suitable for full-time use. This wheelchair has **tilt and recline facility** and is fitted with **elevating leg rests** as standard. This wheelchair has a lightweight Aluminium folding frame and can be disassembled for transportation if required. The seat width and depth can be easily adjusted to perfectly fit the user. This chair offers a wide variety of seating and lying positions for long term use. The **armrests are height adjustable** and can be removed easily as standard and have comfortable polyurethane armrest pads. The **backrest is height adjustable** and is both it and the seat cushion are contoured and padded for excellent support and have a waterproof, breathable and multi-stretch cover.

There is a **solid seat base** which allows use of an alternative cushion if required. The leg rests elevate easily and are fitted with padded calf rests and knee pads; they can also be independently adjusted forward to accommodate uneven leg lengths. **Footplates are height** and angle adjustable as standard. The rear wheels are removable by easy quick release mechanisms and can the chair can be fitted with either Self Propel or Transit wheels. The chair is fitted with **anti-tip wheels** for added safety and a stepper to help an attendant to ease the chair up over kerbs.

Push handles are height adjustable and the **tilt, recline and attendant brakes are clearly labelled** and accessible on the push handles. Mechanical user controlled brakes are also fitted as standard. A padded supportive multi-adjustable headrest and lateral supports as well as clothes guards are all supplied as standard.

- Colour: Metallic Sky Blue
- *See Accessories Section for choice of Armrests & Leg rests Available
- Seat Width: 16" (41cm) to 18" (45cm), 18" (45cm) to 20" (51cm), 20" (51cm) to 22" (56cm)
- Seat Depth: 16" (41cm) to 21" (53cm)
- Seat Height: Adjustable 16" (41cm) - 20" (51cm)
- Backrest Height: 19" (48cm) – 24" (61cm)
- Total Weight: 29kg
- Transport Weight: 16.4kg

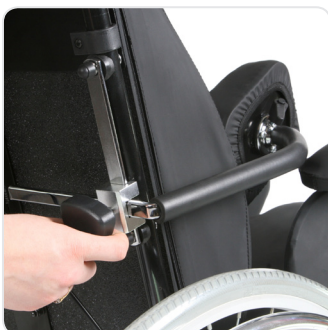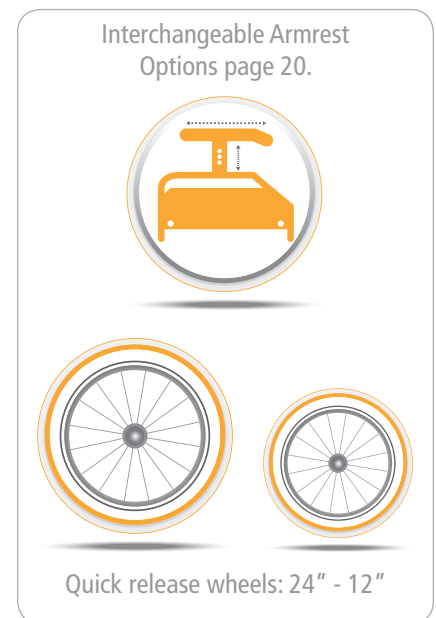

Height adjustable lateral supports as standard.



Width adjustable lateral supports as standard.



Supportive padded seating system as standard.



Pommel removable.
*Optional Accessory



**Max User Weight:
150kg/23.5st.**

X7 Supplied with laterals as standard.



X7 folds for transportation in car boot.



Crash Tested



Puncture Proof Solid Tyres

X8 Series - Paediatric Tilt-in-Space Wheelchair

The **Aktiv X8** is a **paediatric multi-positional** wheelchair suitable for full-time use. This wheelchair is a folding frame and **tilt and recline facility** and is fitted with **elevating leg rests and multi adjustable headrest** as standard. This wheelchair can be even be disassembled for transportation in a car boot if required. The seat width and depth can be easily adjusted to perfectly fit the user. This chair offers a wide variety of seating and lying positions for long term use and a flat position ideal for changing or catheterisation. The **armrests are height adjustable** and flip away as standard and have comfortable polyurethane armrest pads. The backrest is height adjustable and both it and the seat cushion is padded for excellent support and has a waterproof, breathable and multi-stretch cover. A choice of flat or contoured back and seat cushions are available. Both the backrest and seat can be easily removed to allow the chair to fold. The leg rests elevate easily and are fitted with padded calf rests. The footplates are height adjustable and can be fitted with ankle huggers or toe straps if required.

The rear wheels are 16" Transit wheels with **PU puncture proof tyres** and the chair is fitted as standard with anti-tip wheels for added safety and a stepper to help an attendant to ease the chair up over kerbs. **Push handles are height adjustable** and the **tilt, recline and attendant brakes are clearly labelled** and accessible on the push handles. Mechanical user controlled brakes are also fitted as standard. A padded supportive multi-adjustable headrest and lateral supports as well as clothes guards are all supplied as standard. A full range of harnesses and belts available for additional positioning and support.

- Colour: Metallic Sky Blue, Princess Pink , Fire Engine Red, Sunshine Yellow
- Seat Width: 12" - 15" (30- 38cm)
- Seat Depth: 12" - 15" (30- 38cm)
- Seat Height: 17" (43cm)
- Backrest Height: 12" - 15" (30- 38cm)
- Total Weight: 14.8kg
- Transport Weight: 12.4kg

FEATURES

- Tilt function via gas strut
- Reclining function via gas strut
- Folding frame for transportation
- Push handle height adjustable
- Multi-positional headrest
- Removable pommel
- Removable adjustable Lateral supports
- Flip back with padded arm pad
- Swing away elevating leg rests with black plastic foot plate & calf pad
- Aluminium brake
- Breathable upholstery with padded cushion
- 7" grey) PU castor
- Fixed 16" mag rear wheel with grey PU tyre
- Anti-tipper as standard
- Choice of press button or buckle style safety belt





**Max User Weight:
100kg/15.5st.**



Easy transport in a car boot.



Removable pommel.



Seat tilt for maintaining positive postures.



Lays flat for easy changing or resting.



Crash Tested



Puncture Proof Solid Tyres

Product complies with EN12183 and Crash Test ISO7176-19:2001.

Options

Armrest Options



Armrests (suitable for X2, X3, X4, X5, X6)

Our Armrests are interchangeable and flip back and detach as standard. Choose from:



Option 1:
Standard Full Length.



Option 2:
Short Desk Arm.



Option 3:
Loop Arms (grip point
for standing transfers).



Option 4:
Height Adjustable Armrest
(with sliding armpad).

Height adjustable armrest features.

To adjust the height: Simply press the silver button and pull up or push down to change height. Release button to lock position.



Low Height



Higher Position

To adjust the armpad position:

Height adjustable with sliding pad Simply press the button on the side of the armpad and slide forward or back as required.



Desk Length





Full Length

Legrest Release Options (suitable for X3,X4,X5).



Legrest Options (suitable for X3,X4,X5)

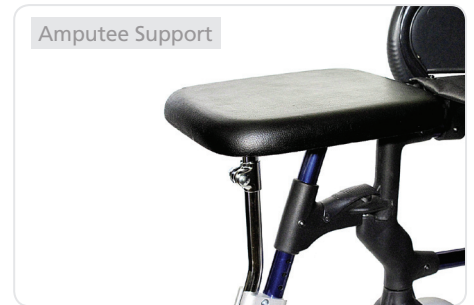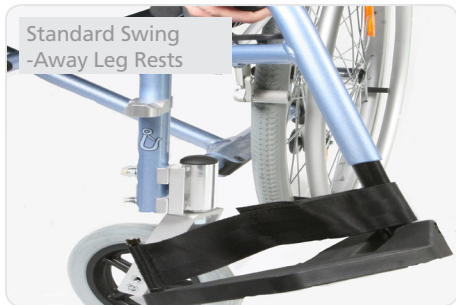

Footplates (suitable for X2,X3,X4,X5)




Anti Tip with Stepper fitted
to x3,x4,x5,x6,x8 as standard.




Options

Universal Options

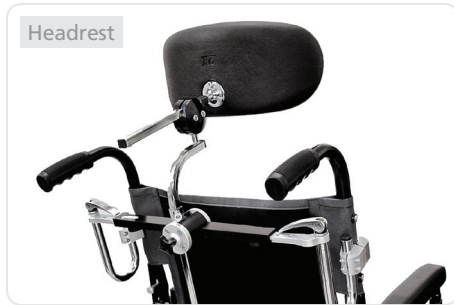

Legrest Options (suitable for X3,X4,X5)

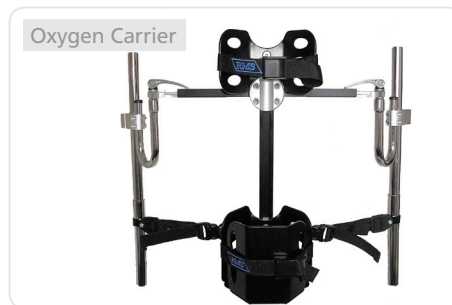

Positioning



WARRANTY

5 Years.

The frame of the wheelchair shall be covered by a 5 year conditional warranty against faults or breakages attributable to faulty manufacture. Definition of the Frame: The frame is the main structure of the wheelchair onto which the various components are fitted to complete the chair. This includes the cross brace and seat base, but not those components, such as the seat that are fitted to it. The frame does not include any removable part of the wheelchair, including those components that require tools for removal. This includes, but is not limited to the entire backrest assembly (including upright canes), Footrests, Armrests, Wheels and Castors, Brakes, Brackets, Upholstery and other Optional Accessories.

1 Year

The various accessories that are fitted to the frame including, but not limited to the entire Backrest Assembly (including upright canes), Footrests, Armrests, Brackets, Upholstery and other Optional Accessories are covered by a 1 year conditional warranty against faults or breakages attributable to faulty manufacture.

CONDITIONS

The warranty commences from the date of purchase. A proof of purchase clearly showing the date of purchase and the dealer's name and details must be produced in order to claim for a warranty repair or replacement. Any warranty claims within the warranty period will be assessed and proven to be defective by an authorised Aktiv Wheelchairs technician before any parts are repaired or replaced. Warranty claims will not be accepted without proof of purchase. This warranty does not extend to freight charges incurred to return the wheelchair to the dealer or to Aktiv Wheelchairs for repair or replacement. The above mentioned warranties are for faults or breakages caused during normal domestic use attributable to only faulty manufacture. This excludes defects or injury caused by or resulting from corrosion, misuse, abuse, neglect, accidental damage, improper assembly, improper voltages, vermin infestation, or any alteration that effects the reliability of the product. Items affected by normal wear and tear are not covered by warranty. This includes, but is not limited to wear on tyres, wheels and bearings, seat and backrest upholstery, armrest pads, push handle grips and brake shoes.

Any alteration that is required to the wheelchair to enhance its function as a request by a qualified Medical Practitioner will require a written endorsement from a Aktiv Wheelchairs engineer. Any alteration carried out without written endorsement will void the warranty. The length of the warranty period may be adjusted at the discretion of a Aktiv Wheelchairs engineer depending on the type of alteration carried out. Any alterations or adjustments carried out without the use of genuine Aktiv Wheelchairs parts will void the warranty. Products will only be replaced under warranty if the fault cannot be repaired. Faulty components that cannot be repaired will not warrant the replacement of the entire product, but will be limited to the replacement of the faulty component only. Every replaced or repaired part shall under the same warranty conditions of the wheelchair.

The warranty period of repaired or replaced parts shall not extend beyond the original warranty period of the wheelchair. Any extended warranties offered by dealers are not part of the warranty offered by Aktiv Wheelchairs. Any claims resulting from an extended warranty arrangement will not be accepted by Aktiv Wheelchairs. The specifications in this manual are general, and may differ slightly from your wheelchair. Warranties cannot be claimed on discrepancies of information in this manual to the specifications of your wheelchair, unless it is proven to be genuine manufacturing fault. Any warranty claims must be processed through the dealer from which the wheelchair was purchased. In the event that a dealer does not satisfy their obligation to process a warranty claim, please contact your agent.



“Mobility, Comfort and Style”



Northern Ireland
Fred Storey Ltd
Old Crow Building
2 Glen Rd - Comber
Co Down - BT23 5EL

TEL: 028 9187 0033
info@fredstorey.com
www.fredstorey.com

MURRAYS

Republic of Ireland
Murray's Medical
Airton Park
Airton Rd - Tallaght
Dublin 24 - Ireland

TEL: 01 8555 696
www.murrays.ie
service@murrays.ie



www.aktiv-wheelchairs.com